


CEC/TH B70A050

Thermostated
 Effective volume: 300 ml to 450 ml
 For disc sample \varnothing : 1.4 cm
 Sample active \varnothing : 1.1 cm
 Sample active area: 0.95 cm²

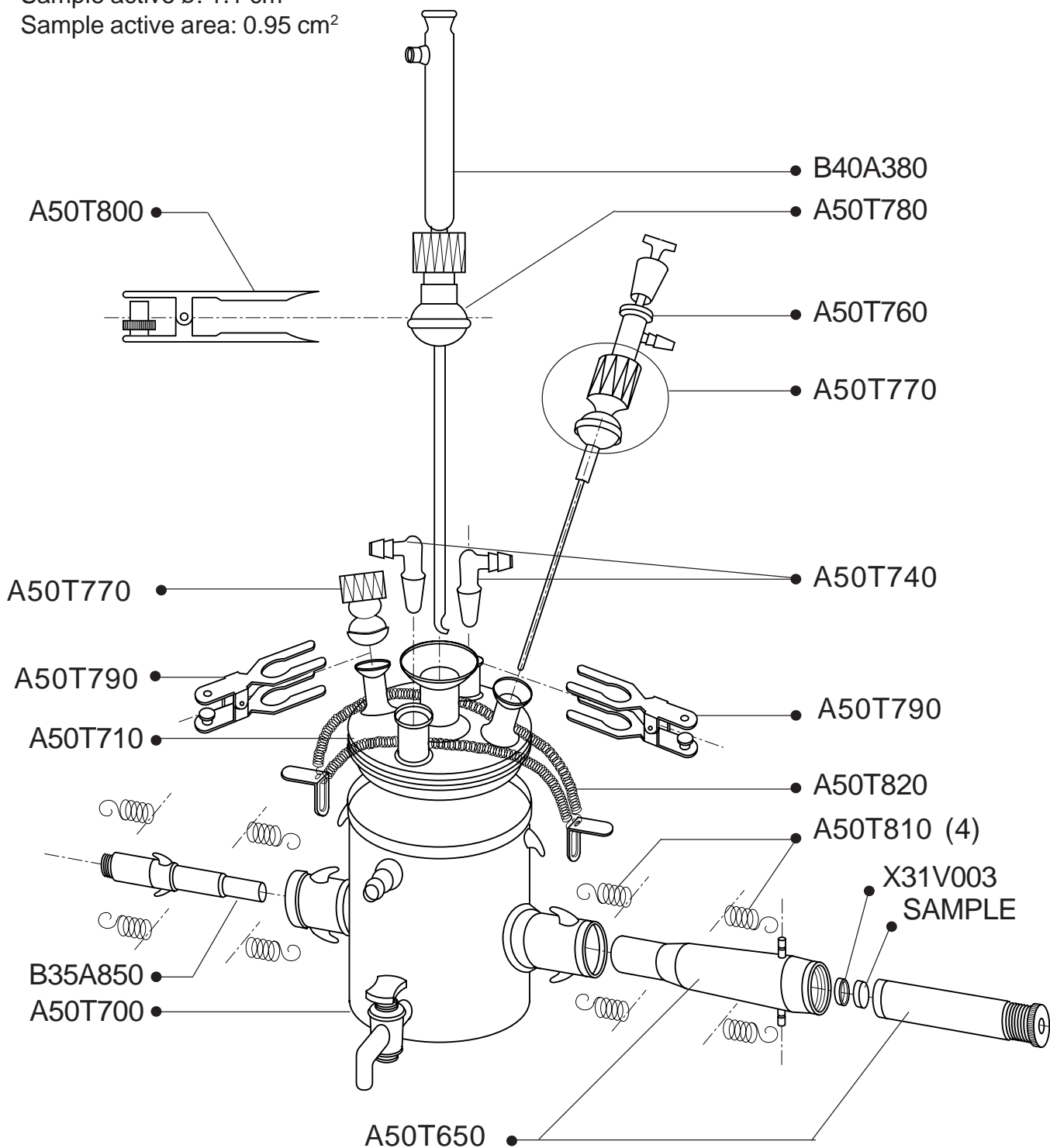

CEC/TH

Thermostated Multipurpose Cell

B70A050

Description

The CEC/TH cell consists of a glass cell with a thermostatic jacket and a disc sample holder offering a reduced active surface making it ideal for small samples. The PCTFE sample holder ensures excellent contact with the sample while avoiding any interaction with the medium. Sample holder and counter electrode are both fitted with ground joints for accurate positioning.

Ordering Information

The CEC/TH product includes

<i>Part No.</i>	<i>Qty</i>	<i>Type</i>	<i>Description</i>
A50T700	1	RMT-RL	Thermostated vessel
A50T710	1	C5R-CS	Cover with 5 ground joints
A50T780	1	RS29/V22/J10	Joint for SVL22
A50T770	1	RS29/V15/J6/J8	Joints for SVL15 (2 pcs)
B35A850	1	CE-PTD-K	Auxiliary Pt electrode for CEC/TH cell
A50T650	1	PEK29-TH	PCTFE sample holder for thermostated cell
B40A380	1	MDL2/CEC	Salt bridge extension for CEC/TH cell without ball and socket joint
A50T740	1	BRT	Tubing for gas bubbler (2 pcs)
A50T760	1	BAROB	Bubbler with stopcock
A50T810	1	RESS/CEC	Springs for CEC/TH cell (4pcs)
A50T800	1	PRS29/15	Clamp for spherical ground joints
A50T790	2	PRS19/9	Clamp for spherical ground joints
A50T820	1	ATT/CEC	Cover holder spring for CEC/TH

CEC/TH consumables

<i>Part No.</i>	<i>Qty</i>	<i>Type</i>	<i>Description</i>
X31V003	1		PCTFE washer for PEK29 sample holder (10 pcs)



RADIOMETER ANALYTICAL SAS
72 rue d'Alsace, 69627 Villeurbanne Cedex, France
E-mail: radiometer@analytical.com Web: www.radiometer-analytical.com
Tel. : +33 (0)4 78 03 38 38 - Fax : +33 (0)4 78 68 88 12